

### DataBank USB Interface and Power Cable (P/N 021-2903)

This cable is included as a standard accessory with the DataBank. This cable is used for communicating and configuring the DataBank as well as charging the internal battery. This cable connects to a computer via USB.



### DataBank RS232 Interface and Power Cable (P/N 021-2900)

This cable is an optional accessory for the DataBank. This cable is used for communicating and configuring the DataBank as well as charging the internal battery. This cable connects to a computer via RS232 Serial port and can be used instead of P/N 021-2903.



### Communication and Integration Cable (P/N 2300-115)

This cable can be used as a standard USB to Serial Adaptor Cable and has been tested using our instruments.



### +12 VDC Power Supply (P/N 7000-941)

This is included as a standard accessory with the DataBank. The power supply connects to the power connector on the interface cable which will supply power to the DataBank. This includes a US line cord P/N 046-0400. If you require another style line cord please specify:

P/N 046-0108 - UK

P/N 046-0125 - Australia

P/N 046-0150 - Europe



### +12 VDC Power Supply Car Adaptor (P/N 2900-151)

This is an optional accessory for the DataBank. This cable connects to the power supply connector on the interface cable and a car cigarette lighter receptacle or other similar receptacle in order to charge the battery of the DataBank.



### External Power Option P/N 2900-940 includes External Power Cable (P/N 2900-152)

This is a factory-installed optional accessory for the DataBank. With the external power option the DataBank can be powered using an external power source such as a 12 volt marine battery. This accessory includes an internal power converter, connector and 3 meter cable assembly. This is recommended for use with the C-Sense



### Cyclops Pigtail Cable with Locking Sleeve


- ◆ 0.6 meter pigtail cable with locking sleeve (P/N 105-0580) shown
- ◆ 5 meter pigtail cable with locking sleeve (P/N 105-0585)
- ◆ 10 meter pigtail cable with locking sleeve (P/N 105-0581)
- ◆ 25 meter pigtail cable with locking sleeve (P/N 105-0582)
- ◆ 50 meter pigtail cable with locking sleeve (P/N 105-0583)



These cables are required for use with the Cyclops-7.

Depending on the P/N ordered, one of these cables comes factory-installed on the DataBank.

Color	DataBank Terminal	Cyclops Sensor
Black	1 GND	2 Power Ground
Blue	2 CLT-A	5 (x10) Gain
Brown	3 CTL-B	6 (x100) Gain
Red	4 PWR	1 Supply Voltage
White	5 IN +	3 Signal Out
Green	6 IN -	4 Analog Ground

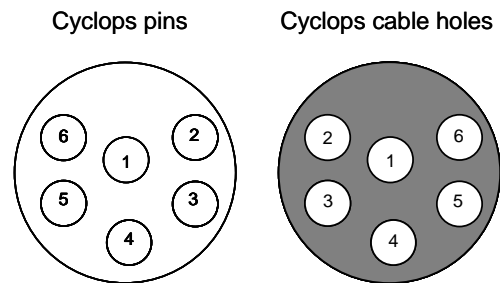

### C-sense Pigtail Cable with Locking Sleeve

- ◆ 0.6 meter pigtail cable with locking sleeve (P/N 105-2410) shown
- ◆ 5 meter pigtail cable with locking sleeve (P/N 105-2411)
- ◆ 10 meter pigtail cable with locking sleeve (P/N 105-2412)
- ◆ 25 meter pigtail cable with locking sleeve (P/N 105-2413)
- ◆ 50 meter pigtail cable with locking sleeve (P/N 105-2414)




These cables are required for use with the C-sense.

Depending on the P/N ordered, one of these cables comes factory-installed on the DataBank.

Color	DataBank Terminal	C-sense Sensor
Black	1 GND	2 Supply Ground
Red	4 PWR	1 Supply Power
White	5 IN +	3 Analog Signal Out
Green	6 IN -	4 Analog Ground

C-sense pins

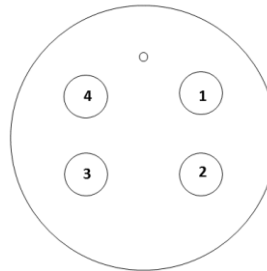

C-sense cable holes

